INSTRUCTIONS FOR USE

Disposable sampling swab

【Product Name】

Disposable sampling swab

【Product Model】

ZCY-03; ZCY-04

【Product structure】

The type ZCY-03, ZCY-04 consist of sampling swab, sheath, catheter, handle and knob. The sampling swab of ZCY-03 consists of a sampling head and a sampling handle, the sampling head is a fluffy head made of PA66 (polyamide 66) material, the sampling handle material is made of ABS (Acrylonitrile Butadiene Styrene) material. The sampling swab of ZCY-04 consists of a sampling head and a sampling handle, the sampling head is made of PU (polyurethane) sponge, the sampling handle is made of PP. The sheath, catheter, handle and knob are made of polypropylene (PP).

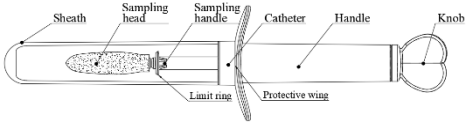


Fig 1 ZCY-03 Structure Diagram

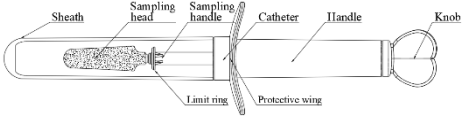


Fig 2 ZCY-04 Structure Diagram

【Intended Use】

It is mainly used for clinical cell sampling and smear.

【Contraindication】

No

【Steps to use】

1. Open the package and take out the sampling swab.
2. Remove the sheath and slowly insert the catheter into the vagina. When the protective wings of the catheter touch the vulvar opening, it is done.
3. Slowly rotate the "heart" -shaped knob clockwise, Push the sampling swab out until the limit ring touches the front end of the catheter.
4. Hold the handle and gently rotate the entire swab clockwise for 3 ～ 5 turns to collect the sample.
5. After sampling is completed, hold the handle and slowly rotate the sampling swab counterclockwise to remove it.
6. After removing the sampling swab, rotate the "heart" -shaped knob counterclockwise to retract the sampling head into the catheter, cover the sheath and send it for testing.

【Precaution】

1. The product is for single use by a single person only. Reuse is prohibited.
2. Please read the instruction manual carefully, pay attention to the expiry date of the product before use. Do not use it after the deadline.
3. Please check the pouch package and the product carefully before use. If either of them is damaged, it is strictly forbidden to use.
4. Before and after sampling, the swab sampling head must not contact other objects to avoid contamination.

5. It is contraindicated for those without sexual history, pregnant women and those during menstruation.
6. 24 hours before sampling, vaginal medication, internal genital irrigation, sexual activity, gynecological examination, and the use of lubricants are strictly prohibited.
7. If you are undergoing treatment for cervical diseases, you must use it under the guidance of a doctor.
8. If obvious pain or bleeding occurs during the sampling process, stop using it immediately.
9. Waste disposal of disposable sampling swab should be in accordance with the sanitary management regulations of the medical.

【Definitions of Signs】

	Sterilized using ethylene oxide		Do not re-use
	Do not use if package is damaged		Do not resterilize
	Caution		Fragile, handle with care
	Upward		Keep dry
	Batch code		Use-by date
	Manufacturer		Authorized Representative in the European Union
	Date of manufacture		Consult instructions for use
	Temperature limit		Keep away from sunlight

	Stacking limit by number	MD	Medical device
CE	" CE " Marking	UDI	Unique Device Identification

【Storage and transportation】

1. The product should be stored in a room with a temperature of -10℃～40℃, relative humidity of not more than 80%, good ventilation, no corrosive gas.
2. The product should be protected from heavy pressure, strong vibration, direct sunlight, high temperature, moisture, rain and corrosive gases or liquids during transportation.

【Sterilization】 EO sterilization

【Shelf Life】 5 years after sterilization



Name: Jiangsu HanHeng Medical Technology Co., Ltd.
Add: #118 Xinyuan 3rd Road, Xinbei District, Changzhou, 213031 Jiangsu P.R. China
Tel: 0086-519-85559288
Web: https://hanheng-medical.com

EU REP

Name: Luxus Lebenswelt GmbH
Add: Kochstr. 1, 47877, Willich, Germany
Contact Person: Liu Sun
Tel/Fax: 0049-1715605732
E-mail: Info.m@luxuslw.de

【Version】 2025.03.26 A/0

CE 0197